


danuTips for Enteral Care



How Polyurethane and TRITAN™ are Advancing Enteral Devices

In the medical devices industry, excellence and innovation are paramount. We prioritize using top-tier materials to ensure product longevity and safe application.

Polyurethane (PU) →

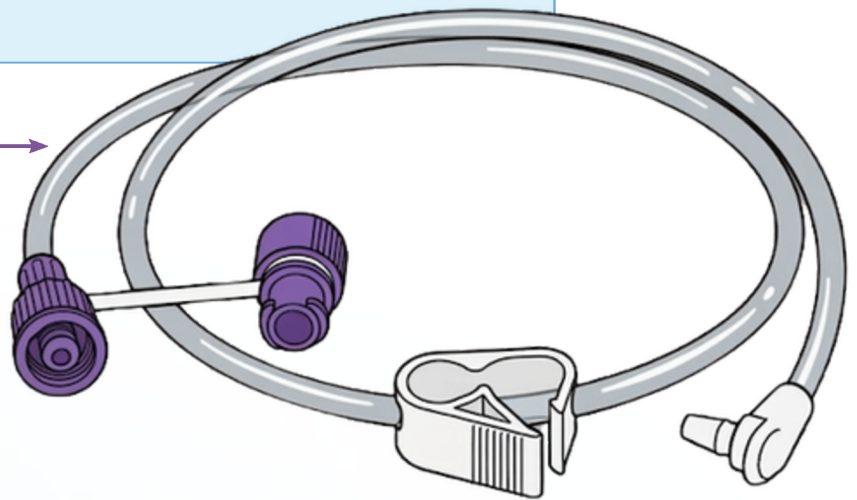

Polyurethane (PU)


We choose Polyurethane over PVC for the tubing in our danumed® Extension Sets due to its exceptional qualities:

- ✓ **High biocompatibility:** PU is a biocompatible material, reducing the risk of adverse reactions and irritation, and making it ideal for long-term use.
- ✓ **No need for plasticizers:** Unlike PVC, PU remains flexible and soft without harmful plasticizers like phthalates, preventing potential health risks for patients.
- ✓ **Resistance to kinking:** PU tubes are less likely to kink which is crucial for maintaining a consistent flow in enteral feeding.
- ✓ **Durability and flexibility:** PU retains its flexibility and elasticity over time while offering excellent tensile strength and elongation.
- ✓ **Enhanced flow range:** PU tubes have thinner walls than PVC tubes, allowing for a larger inner diameter and enhanced flow area for enteral feeds.
- ✓ **Patient comfort:** PU's softness and flexibility provide superior comfort, particularly for long-term or sensitive applications.

TRITAN™

TRITAN™ is our material of choice for the core pieces of our enteral connectors because:

- ✓ **Durability:** It's a lightweight, robust plastic known for its resistance to aggressive solutions and medicines.
- ✓ **Stress crack resistance:** TRITAN™ reduces the risk of stress cracks, supporting break-resistant use.
- ✓ **Structural integrity:** It can withstand the rigors of daily use without compromising the device's structural integrity.

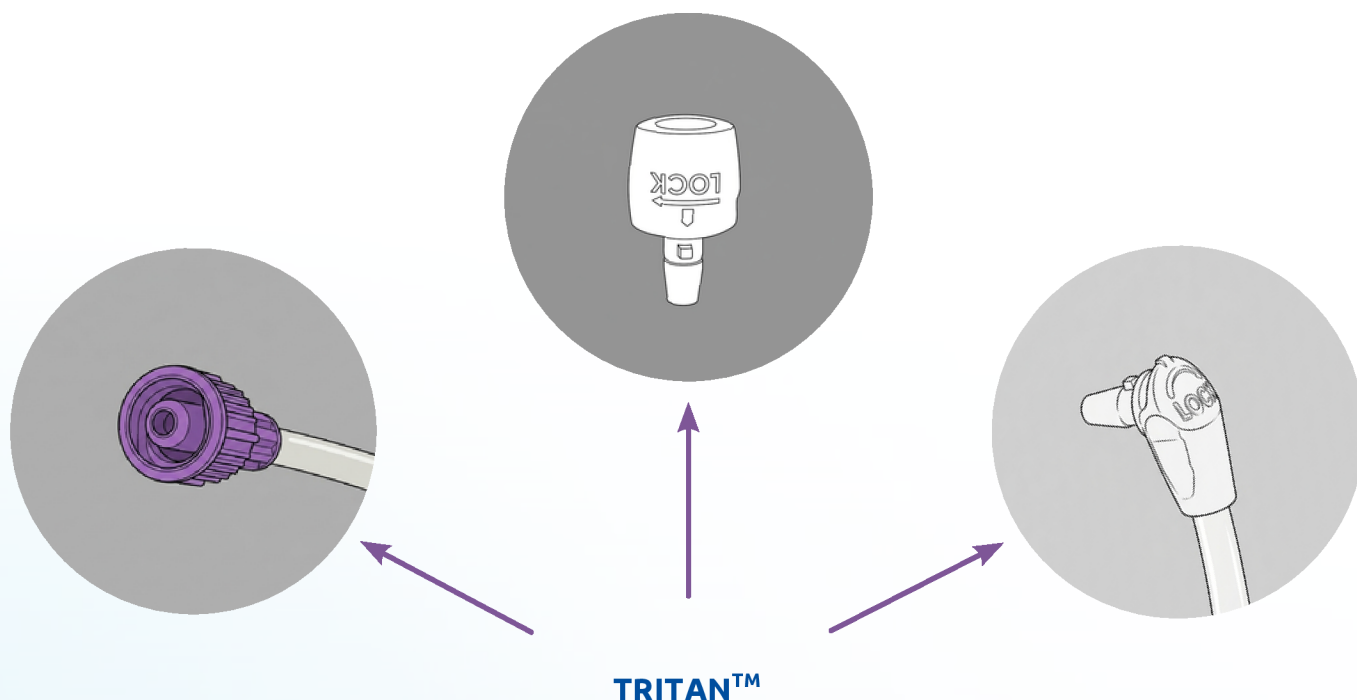

By carefully selecting materials like TRITAN™ and Polyurethane, we enhance the reliability, safety, and longevity of our danumed® products. These choices are crucial to protecting the health of vulnerable patients who depend on safe enteral applications for their nutritional and medical needs.

Disclaimer: This information is not a substitute for professional medical care. In case of problems, complications, or questions, always contact your medical professional.

References:

Dual-purpose gastric decompression and enteral feeding tubes rationale and design of novel nasogastric and nasojejunal tubes; Silk David, Quinn David; *Journal of Parenteral and Enteral Nutrition*, Vol. 39 No. 5, 2015
 Biomedical Applications of Polyurethanes: A review of past promises, present realities and a vibrant future; Zdrahala Richard et. al., *Journal of Biomaterials Applications*, Vol. 14, 1999
 Biomedical Applications of Polyurethanes; Bergeron Mylène, Guidoin Robert et. al., Chapter 8, 2001
 Carolyn Best (Nutrition Nurse Specialist, Royal Hampshire County Hospital, Winchester); *Nursing Times (online)*, March 2019, Vol 115 Issue 3
 Ösophagusperforation und Pneumothorax - Komplikationen durch Legen einer Magensonde (Fallbericht); Tronnier V. et. al., *Anästhesiol. Intensivmed. Notfallmed. Schmerztherapie* 26, 1991
 Migration von Weichmachern aus PVC Schläuchen in enterale Nahrungslösungen; Welle Frank et. al., *Pharma International* 3/2005
 Evaluation of the elastic behaviour of central venous PVC, polyurethane and silicone catheters; Cervera M et. al., *Phys. Med. Biol.* Vol. 34, 1989
 Chemical compatibility with hospital disinfectants and oncology drugs, Eastman Brochure, September 2014
 Lipid and Isopropanol Resistance of Eastman Polymers in Medical Devices, Eastman Brochure, September 2007
 Redefining the balance between processability and chemical resistance, Eastman Brochure, February 2009
 Illustrations created with AI assistance

SCAN ME FOR
MORE TIPS

